


# LIS Basics: CP and AP LIS Design and Operations

**Pathology Informatics 2014**  
**May 13, 2014**  
**Pittsburgh, PA**

Walter H. Henricks, M.D.  
Cleveland Clinic

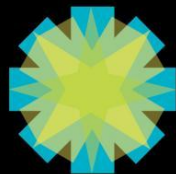

**P A T H O L O G Y**  
**I N F O R M A T I C S**  
**S U M M I T 2 0 1 4**

May 13-16, 2014  
Pittsburgh, PA


Brought to you by the Association for Pathology Informatics.

# Introduction to Laboratory Information Systems (LISs)

## Learning Objectives

- Define terms and jargon related to LISs
- Describe central importance of dictionaries and interfaces to LIS function and laboratory operations
- Identify LIS functions as they relate to laboratory workflow in CP and AP

# The LIS is the Engine of the Laboratory

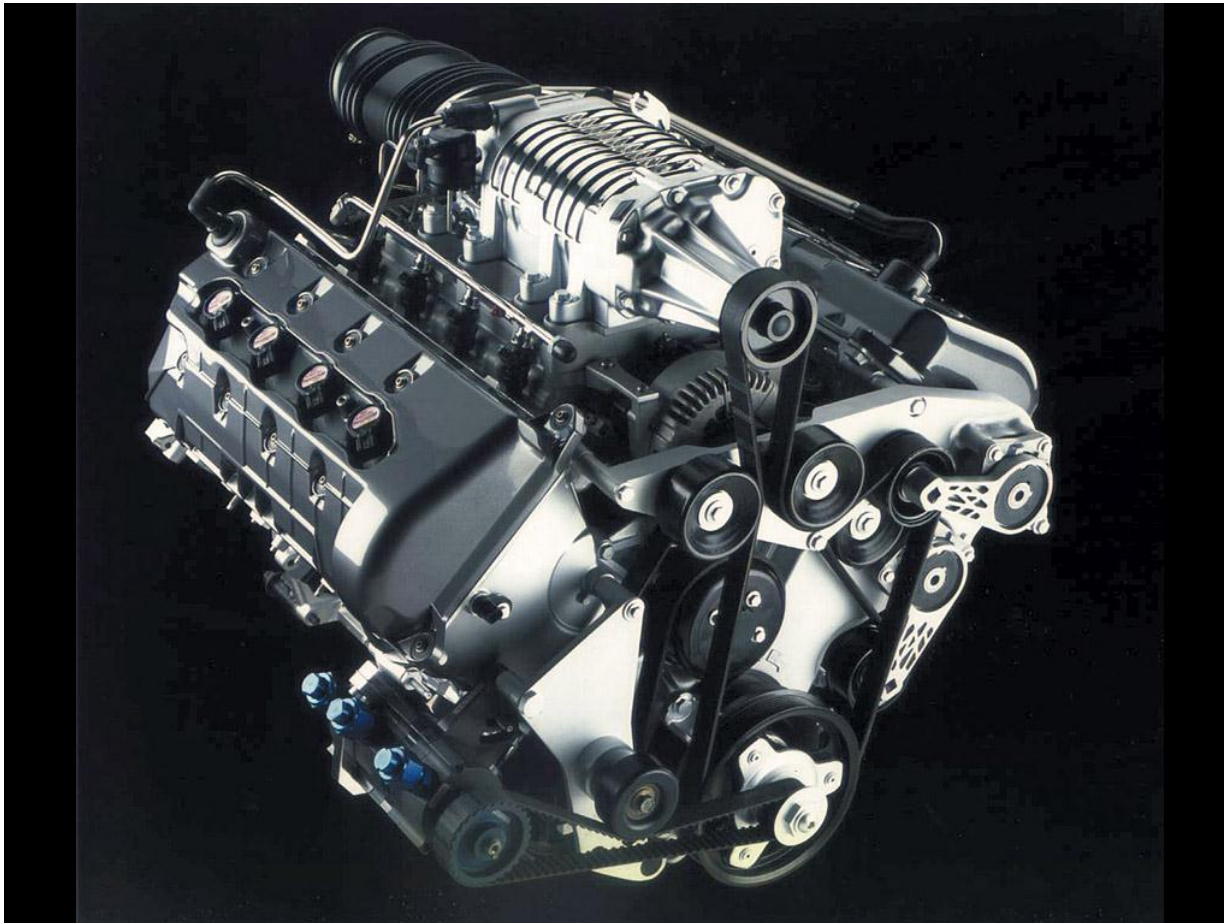

# Laboratory Information System (LIS)


- Interrelated programs and hardware that provide electronic data processing and information management functions necessary for laboratory operations
- Database that establishes and maintains standard definitions and information processing procedures (Elevitch and Aller)

# Introduction to LISs – Outline

## ➤ LIS architecture

- LIS dictionaries (a.k.a. maintenance tables)
- LIS functions in laboratory workflow
  - Clinical laboratory (CP)
  - Anatomic pathology (AP)

# Building Blocks of Laboratory Information Systems

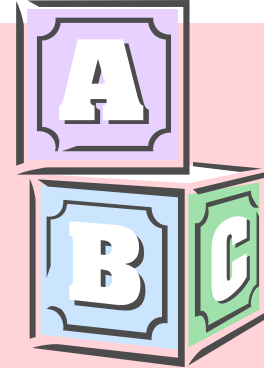

LIS application software

Database Management System  
(DBMS)

Operating System


Hardware





Peter B. Lewis Building, Case Western Reserve University. Architect: Frank Gehry

# Host-based LIS Architecture

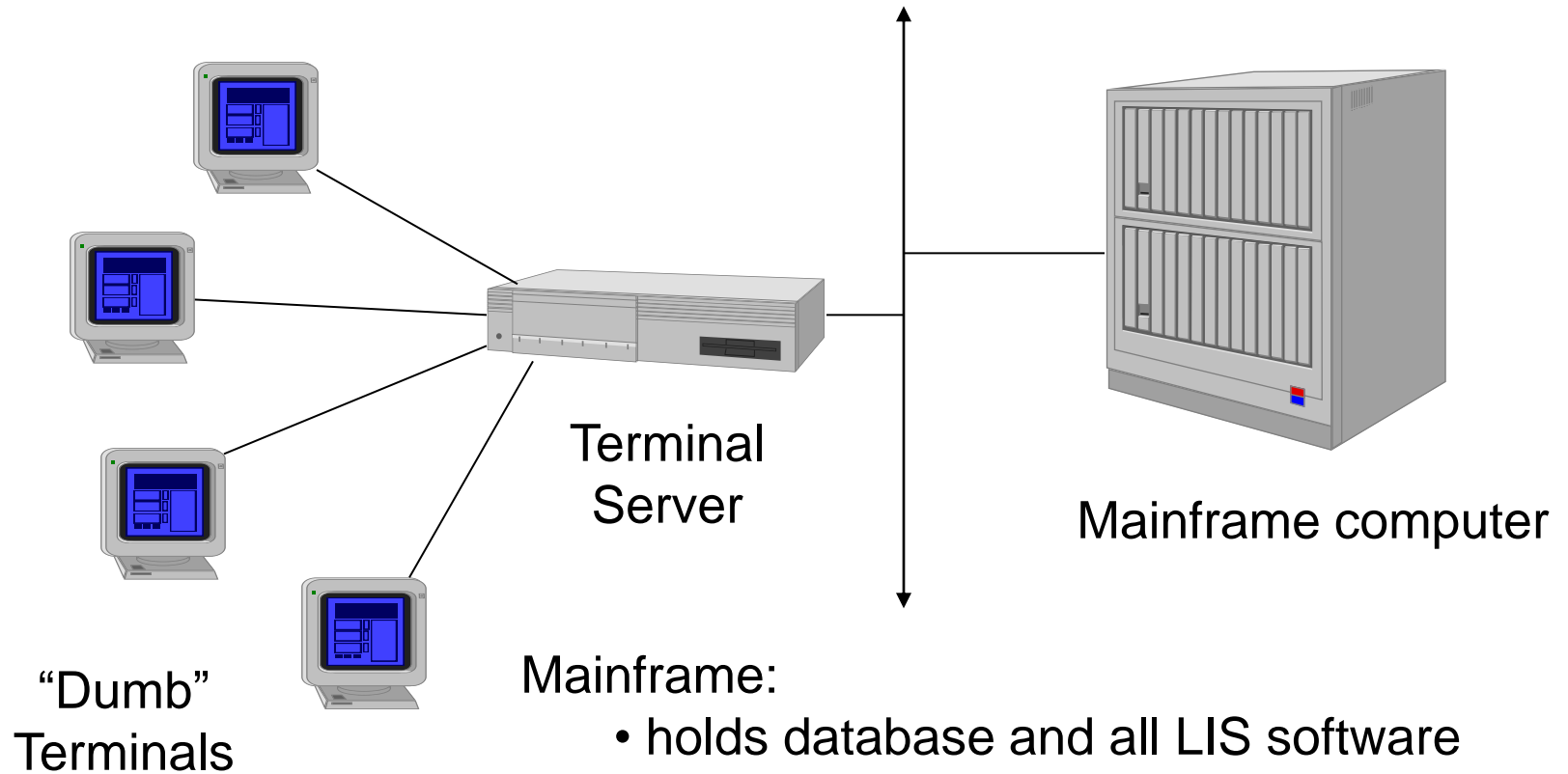

Mainframe:


- holds database and all LIS software
- manages all LIS functions and transactions

Terminals:

- data display and input only
- PCs can connect using "terminal emulation"



# Client/Server LIS

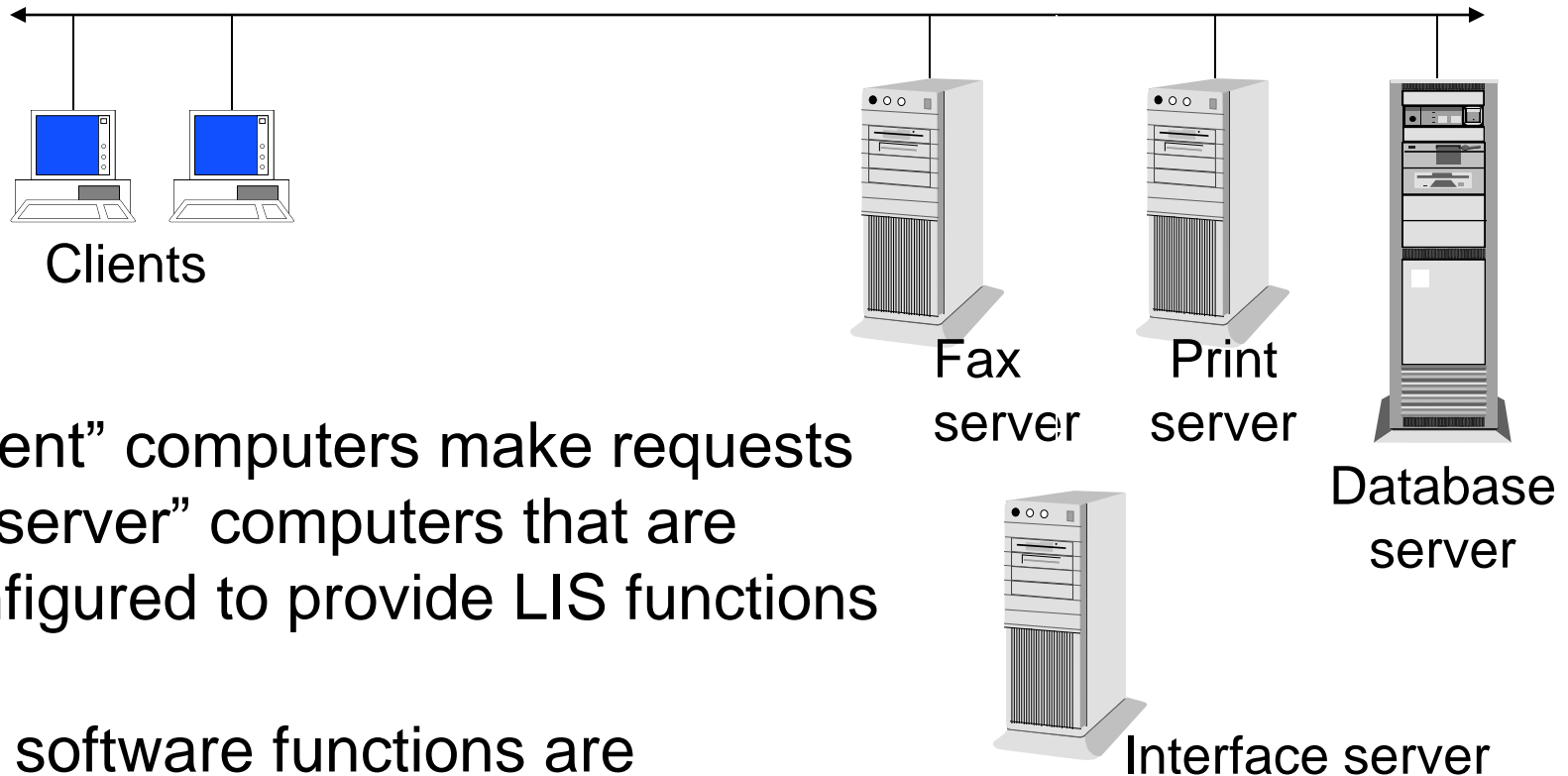

“Client” computers make requests of “server” computers that are configured to provide LIS functions

LIS software functions are **distributed** across all clients and servers

# Mainframe vs. Client/Server LIS




- **Mainframe/Host-based**

- Typically character-based user interface
- Limited flexibility
- Good security
- Centralized maintenance and control
- Single point of failure

- **Client/Server**

- Graphical user interface
- Greater configurability
- Greater security risks (e.g. viruses, PC ports)
- Decentralized and more distributed maintenance
- Multiple points of failure, though each less catastrophic

# Thin Client Architecture in LIS

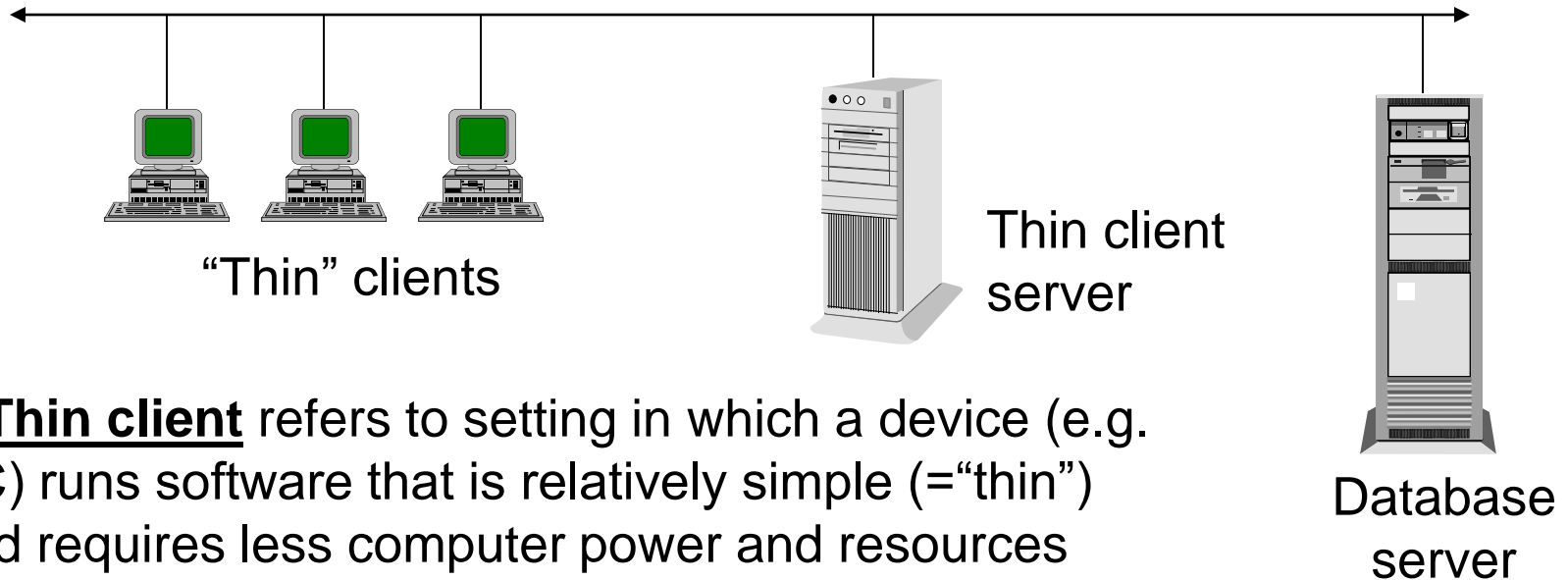

- **Thin client** refers to setting in which a device (e.g. PC) runs software that is relatively simple (=“thin”) and requires less computer power and resources

- Application logic executes on *thin client server*
- Resembles host-based/mainframe model in some respects (connection through intermediate server)



# Thin Client Computing for the LIS

## How it may benefit **YOUR** laboratory

- Easier administration
  - Standardized application/programs controlled centrally
  - Easier to implement software updates in a complex environment
- Cross-platform (PC, Mac)
- Lower hardware requirements and costs
- Remote access
- Less network traffic

# Thin Client Computing for the LIS

## Why **YOUR** laboratory may think twice




- Hardware and license costs
- Single point of failure for all workstations connected to thin client server
- Effectiveness of vendor's implementation of thin client
- Potential inability to do specialized functions on thin client workstation, e.g. imaging, voice recognition

# Introduction to LISs – Outline

- LIS architecture
- LIS dictionaries (a.k.a. maintenance tables)
- LIS functions in laboratory workflow
  - Clinical laboratory (CP)
  - Anatomic pathology (AP)



# What do dictionaries contain?

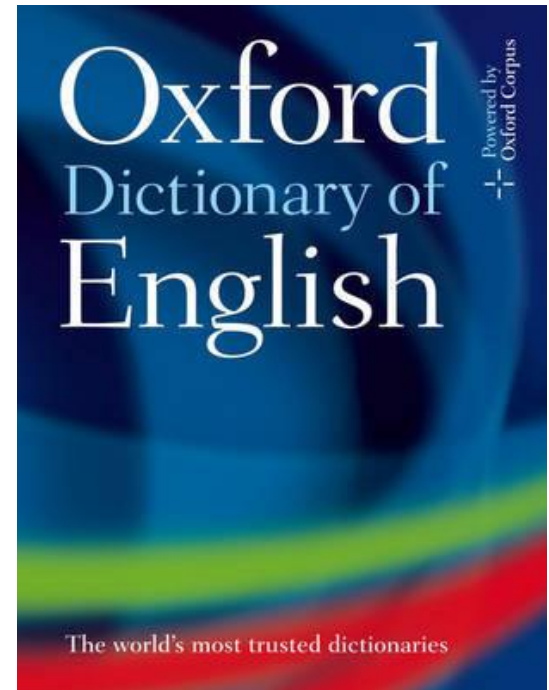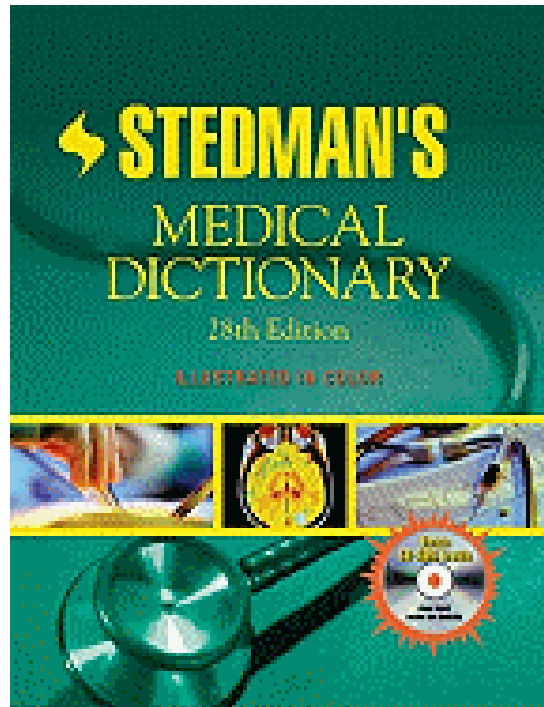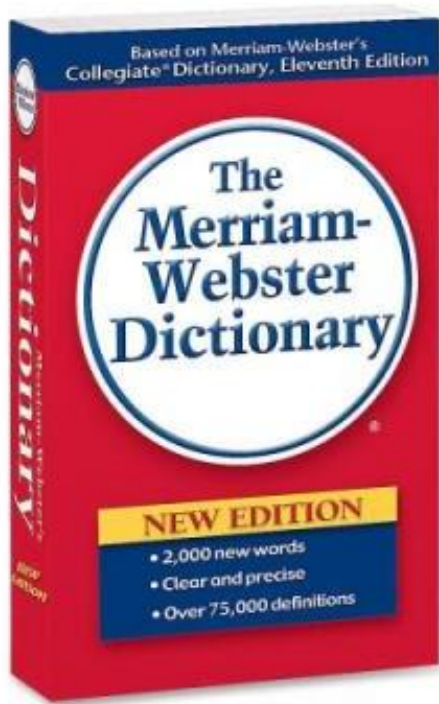

## DEFINITIONS!

# LIS Dictionaries *Define* Your Laboratory's Information Framework

- Standardize naming conventions and procedures
- Standardize laboratory terminology and definitions
- Constrain choices for data fields to ensure valid data entry
- Define content and format of reports (e.g. units of measure)
- Define rules and calculations

# LIS Dictionaries Define *People, Places, Things*

- Test and test battery/profile definitions
- Test worksheets / worklists
- Person dictionaries (e.g., ordering physician, pathologist, technologist)
- Security/access level privileges for user types
- Patient locations
- Laboratory locations (e.g. sections, “areas”)
- Specimen types
- Histologic stain protocols (e.g. Giemsa on gastric bx)
- Analyzer/instrument interfaces
- Autoverification parameters
- Many others...

## TEST DEFINITION DICTIONARY

TEST NAME: Hemoglobin

TEST CODE: HGB

LAB. DEPT: CORE

CONTAINER TYPES: LAV

WORKSHEET(S): CELCOUNTR

IN BATTERIES: CBC, CBCDIF, HGBHCT

AUTOVERIFY RANGE: 6.1-19.9

## LAB DEPT DICTIONARY

CORE  
CLINIC  
GASLAB  
Etc.

## CONTAINER TYPE DICTIONARY

LAV  
BLUE  
RED  
Etc.

## TEST BATTERY DICTIONARY

CBC  
CBCDIF  
HGBHCT  
PTINR  
BMP  
Etc.

## AUTOVERIFICATION DICTIONARY

RULES FOR HGB  
Etc.

## INSTRUMENT INTERFACE TABLE

BLDCTR INTERFACE  
MAINTENANCE  
CHEM INTERFACE  
MAINTENANCE  
Etc.

# Example AP LIS Dictionary: Specimen Part Type

PART TYPE DICTIONARY (mock)	
Entry Name:	LUNG, TXP BX
Shorthand:	TLBX
Description:	LUNG TRANSPLANT BIOPSY
Synonyms:	LUNG
	LUNG BIOPSY
	PULM
Specimen Categories:	SURGICAL ROUTINE
	SURGICAL OUTSIDE REVIEW
	SURGICAL CONSULTATION
Protocol:	LUNG TRANSPLANT BX
Fee Code(s):	LEVEL V

PROTOCOL DICT.
LUNG BX
LUNG TRANSPLANT BX
PROSTATE BX
⋮

FEE CODE DICT.
LEVEL I
⋮
LEVEL VI

## PART TYPE DICTIONARY (mock)

Entry Name:	LUNG, TXP BX
Shorthand:	TLBX
Description:	LUNG TRANSPLANT BIOPSY
Synonyms:	LUNG
	LUNG BIOPSY
	PULM
Specimen Categories:	SURGICAL ROUTINE
	SURGICAL OUTSIDE REVIEW
	SURGICAL CONSULTATION
Protocol:	LUNG TRANSPLANT BX
Fee Code(s):	LEVEL V

## HISTO PROTOCOL DICTIONARY (mock)

Entry Name:	LUNG TRANSPLANT BX
Shorthand:	TLBX
:	
Components:	H&E, INITIAL
	H&E, LEVEL
	MOVAT STAIN
	GMS STAIN

## STAIN DICTIONARY (mock)

Entry Name:	MOVAT STAIN
Shorthand:	MOVAT
:	
Label Print:	MOVAT
Stain Fee Code(s):	SPECIAL STAIN GRP 2


# Dictionaries and maintenance tables *tailor* the LIS to **YOUR** laboratory

- Table definition is critical to successful LIS implementation and requires planning and user involvement.
- LIS “out of the box” may have default entries, which will most likely need to be re-defined.
- Definition changes require careful attention to prevent unintended consequences.
- Changes must be tested before being “put into production”.





# Interfaces are important to life...and to laboratories

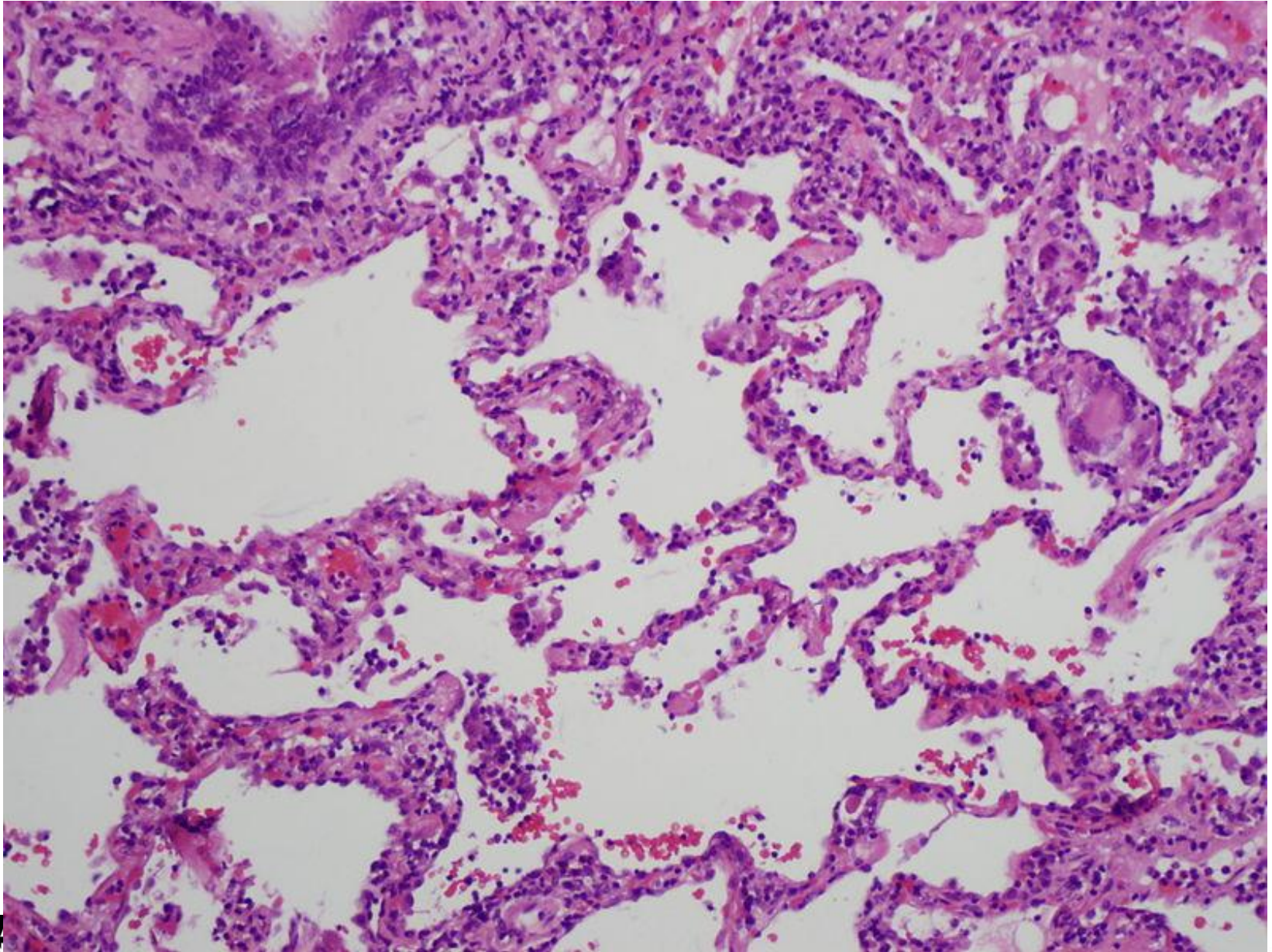

Walter H. Henricks, M.D.



# LIS Interfaces

**Interface** – software and connections that translate electronic messages so that otherwise incompatible systems can *exchange* data

**LIS interfaces are critical to laboratory success (e.g. test order receipt, results reporting)**

# LIS-Instrument Interfaces

- Download = direct transfer of patient identification and test order data from LIS to instrument
- Upload = direct transfer of results back to LIS
- Uni-directional vs. Bi-directional
- Broadcast vs. query
- Unique specimen number on *LIS-generated barcode specimen label* links order and result data in the analyzer and the LIS.
- Other LIS interfaces: handheld phlebotomy devices, tissue cassette engravers, point of care testing devices



# LIS-Instrument Interface Implementation

- LIS vendors have “off-the-shelf” interfaces for most common instruments (revenue source).
- Installation of a new interface is not “plug and play.”
- Interface software must be installed in LIS dictionaries.
  - Definition of data and sequence in the manner expected in the relevant worksheet(s).
  - Rigorous testing and validation prior to use and with changes.

# **YOUR LIS is interfaced to many other computer systems**

- Electronic health record (EHR)
- Admission-Discharge-Transfer (ADT) – pt. registration
- Web portal system
- Physician office systems
- Billing
- Other LISs
- Others...



# HL7 (Health Level 7) – Most Important Data Exchange Standard In Healthcare

- *HL7 defines the format (syntax, structure) but not the specific content of messages*
- HL7 tells computer systems “how to say it” to each other but not “what to say”.

# HL7 Basics

- HL7 defines various message types, such as laboratory test orders and results, patient admission-discharge-transfer (ADT), others.
- Messages (e.g. order, result) consist of segments comprised of fields in a defined sequence and format

```
OBX|1|NM|CK^CK|1|251  
|U/L|30-220|H|||F|||||201306082152|
```





**ZÁKAZ TLUMOČENÍ**  
**TRANSLATING**  
**PROHIBITED**

# HL7 does not eliminate the difficulties of implementing interfaces

- HL7 interface specifications for LIS typically do not match other vendors/systems.
- Institution-specific HL7 segments often exist.
- Translation tables are necessary to cross-reference different test codes in different systems.
- **LOINC** is a standardized vocabulary of lab test names with goal of interoperability
- Interface deployment requires testing, validation, documentation

# Introduction to LISs – Outline

- LIS architecture
- LIS dictionaries (a.k.a. maintenance tables)
- LIS functions in laboratory workflow
  - Clinical laboratory (CP)
  - Anatomic pathology (AP)

# Preamalytic Phase Information Management

- Order creation and test selection
- Specimen collection and labeling
- Specimen receipt and tracking

# Order Creation and Test Selection – CP LIS

- Orders interfaced from EHR
  - Laboratory test orders are entered in EHR (based on menu choices and EHR definitions).
  - Orders cross HL7 interface to LIS.
  - Interface design ensures that correct orders are filed in LIS based on matching or translation of test codes.
- Non-interfaced test requisition order entry
  - Paper accompanies specimen to laboratory.
  - Lab staff enters orders into LIS, with test choices based on LIS dictionaries.

# Specimen Collection and Labeling – CP LIS

- Inpatient – phlebotomist sweep collection
  - LIS assigns unique (accession) specimen and/or container number to each order.
  - LIS places orders on a collection list for next scheduled sweep or “AM labs”.
  - Phlebotomist may generate label dynamically at bedside using portable device.
  - Collection list or labels may guide phlebotomist as to appropriate container type to use.
  - Phlebotomist applies LIS-generated label to specimen.

# Specimen Collection and Labeling – CP LIS

Outpatient – phlebotomy draw station (service ctr.):

- EHR-interfaced orders – phlebotomist/tech accesses existing orders in LIS, prints labels, and collects specimen.
  - There may be a preliminary step that involves locating patient orders in EHR and “pulling them into” LIS
- Non-interfaced orders – phlebotomist/tech enters orders for specified tests into LIS, prints labels, and collects specimen.



# Specimen Receipt and Tracking – CP LIS

- EHR-interfaced orders
  - When specimen arrives in lab, orders already exist in LIS.
  - Lab staff acknowledges specimen receipt in LIS and confirms label and orders.
- Non-interfaced orders
  - Orders do not exist in LIS when specimen arrives.
  - Lab staff orders tests in LIS, prints and applies labels.
- Specimen status in LIS: “received” or “in-lab”

# Specimen Receipt and Order Creation – AP LIS

- LIS assigns specimen a unique accession number (“accessioning”).
- Different “number wheels” can be used to distinguish different classes of specimens
  - e.g. by location, HS-11-123, CS-11-123
- Multiple specimen parts are identified under one accession number.
- Other data from requisition are entered– e.g. physician, clinical hx
- Patient demographic data may be pre-populated in LIS from ADT interface data.

# Specimen Receipt and Order Creation – AP LIS

- Two field types in LIS identify each part:
  - Part type – selected from a dictionary; categories of specimen types, e.g. colon, polyp
  - Part description – free text additional description provided with specimen, e.g. “large colon polyp at 50 cm”
- Part types may be linked in LIS to histology protocols or special stains.
- LIS may print specimen label with bar code.
- Specimen status is updated, e.g. to “Accessioned”.

# Analytic Phase Information Management

- Work distribution and specimen preparation
- Test performance and analysis
- Test interpretation
- Additional testing based on initial results
- Results entry

# Work Distribution – CP LIS

- Many “orderable” tests in CP consist of groups (or panels) of multiple individual test components, e.g. BMP, CMP
- LIS files orders for individual test components.
- LIS-generates bar code labels with unique specimen IDs that are key to instrument interfaces and specimen tracking.
- Test orders are *routed* to the appropriate LIS worksheets based on the worksheets assigned in the LIS test definitions.

# Work Distribution and Worksheets – CP LIS

- Tests performed on interfaced instruments have LIS worksheets linked to instrument maintenance dictionaries, ensuring download to appropriate instruments
- Download to instrument may occur based on different triggers:
  - Receipt of order in LIS from EHR interface
  - Receipt of specimen in lab (as tracked in LIS)
- For batched tests, orders are routed to the appropriate LIS worksheets, which technologists access (or print) to see the list of work for that run.

# Test Performance and Result Entry – CP LIS

- For interfaced instruments:
  - Instrument software reads bar code specimen number and performs the tests per downloaded orders (from LIS) linked to that specimen number.
  - Results are uploaded back to LIS, tied to specimen number
  - Interface specifications shared between LIS and instrument software ensure data transfer in expected sequence and format.
  - LIS worksheets linked to the instrument maintenance ensure that test results are filed correctly in the LIS.

# Test Performance and Result Entry – CP LIS

- For tests performed on non-interfaced analyzers or manually, technologists enter results into LIS worksheets using the LIS resulting function.
- *Footnotes or comments* may be required to add additional information – free text vs. coded template from LIS dictionary



# Rules and Additional Testing – CP LIS

- Worksheets link test or battery to any rules or calculations to be performed based on initial results, e.g. anion gap.
- Reflex Testing – automatic generation of new test order in LIS based on initial results meeting defined criteria, e.g. titration of positive ANA screen
- Autoverification – automatic final verification in the LIS of results from automated instruments without manual intervention.
- Criteria are based on algorithms defined in LIS
- Results or specimen-related data from the instrument that fall outside defined criteria are held for review.

# Autoverification Table in LIS

TEST: HGB

Use Normal Range (<Y>/N)	: N	
Use Borderline Range (<Y>/N)	: N	
Use Technical Range (<Y>/N)	: N	
Use Verify Range (<Y>/N)	: N	
Use Delta Check (<Y>/N)	: N	
Use Instrument Filing Range (Y/<N>)	: Y	(Fail Cup)
Use Invalid (???) Range (Y/<N>)	: N	
Fail on Result Flag(s) (Y/<N>)	: Y	Include: 4,
Fail on Pattern(s) (Y/<N>)	: N	

# Middleware

- **Middleware**: rules-based processing provided by instrument vendor or third party that “sits between” the LIS and instrument
  - Autoverification
  - Reflexive test ordering based on result
  - Automatic dilutions, repeats, smear creation
  - Other aspects of instrument management, e.g. maintenance alerts

# Grossing and Specimen Preparation – AP LIS

- Main information outputs in LIS of “grossing” phase are:
  - Text entry in LIS “Gross Description” field
  - Use of pre-defined templates for some specimen types
  - Designation of tissue sections in Histology module
  - Status updated, e.g. “Gross Complete”

# Slide Preparation and Work Distribution – AP LIS

- LIS directs slide preparation workflow
  - Pre-defined protocols for levels and stains
  - Histology logs defining worklists of cases and blocks from grossing step
  - Slide labels (+/- bar codes) based on data entered in histology module and protocol/stain definitions
  - Special stains appear on specified logs (e.g. immunohistochemistry log)
- “Asset” tracking based on bar code labels and points of scanning defined

# Slide Preparation and Work Distribution – AP LIS

- LIS produces “working draft” report for pathologist – paper vs. paperless
- Working draft format and content are based on template configuration in LIS, e.g.
  - Clinical information
  - Gross description
  - Frozen section report (if performed)
  - Summary of patient’s previous results, based on LIS search of database

# Report Generation – AP LIS

- Entry of Final Diagnosis in LIS facilitated by:
  - Pre-defined templates, checklists, and formats
  - Speech to text conversion capability
  - Automatic entry of billing (CPT) and diagnosis (ICD) codes based on dictionary definitions
- “Synoptic” LIS modules enable entry of *discrete data elements*.
- Pathologist signs out cases with electronic signature that locks the case.

# Post Analytic Phase Information Management

- Generation and delivery of lab results and test reports
- Correcting, amending, and addending reports



# Report Distribution

- Hard copy report format is based on configurable template in LIS.
  - Printing – scheduled batches, on demand
  - Faxing – automated if fax numbers in dictionary
- For electronic reporting, reports pass from the LIS to receiving system via an interface
- *The format and display of interfaced reports is dictated by the screen design in the receiving system.*
- PDF and RTF interfaces preserve formatting; receiving system must accommodate.

# Corrected Results – CP LIS

- LIS report must clearly identify the corrected result as corrected.
- Corrected result must also include the original result.
- Corrected result typically also includes documentation of the person correcting the result and a record of any communications (e.g. “corrected result called to ...”).
- Corrected report typically replaces (overlays) previous result in EHR; original kept in audit trail.

# Amendments and Addenda – AP LIS

- Report formats should be configured so that amended or addended status is obvious.
- Entire report is re-printed or re-transmitted across the interface with the new addendum identified as such
- In the EHR, the new report overlays the previous report.

# Microbiology – LIS Considerations

- Specimen sources (types)
- Specimen description and clinical data
- Organism dictionaries
- Preliminary reports, sequential updates
- Antibiotic susceptibility testing and reporting

# Blood Bank/Transfusion Medicine LIS Considerations

- Blood component identification
- Component status – e.g. allocated, issued
- Antibody work ups
- Special needs – e.g. irradiation, CMV neg, etc.
- Donor center functions
- Transfusion reactions
- Apheresis
- FDA approval; validation; inspections

# LIS Fundamentals – Summary

- LIS dictionaries define the framework for information processing and workflow.
- Worksheets and logs define the data and specimen flow in the laboratory.
- LIS is central to laboratory operations throughout all phases of testing.
- Capabilities in LISs reflect workflow differences between CP and AP.